




StealthStation[®] S7[®] System



StealthStation® S7® System

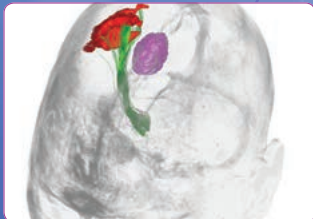

RIGHT INFO, RIGHT TIME

Information is vital in the Operating Room. And the StealthStation® S7® System is built to put the intraoperative imaging and surgical navigation information you need in the right place, at the right time.

- **The Bigger Picture.** The StealthStation® S7® system comes with a 24 inch (61 cm) high resolution widescreen display, as well as the option to make that a touchscreen display.
- **Intraoperative Imaging's Partner.** It interfaces with Medtronic's versatile 3D imaging and 2D fluoroscopy O-arm® system and Polestar® surgical MRI system, as well as iMRI, iCT, and C-arm imaging systems.
- **State-of-the-Art Tracking Options.** The new compact design houses both Medtronic tracking technologies: the advanced optical navigation camera and our unique AxiEM™ electromagnetic system – giving you the choice to navigate each procedure the way you want.
- **Easier Access to Pre-Op Exams.** Easier access to Pre-Op Exams with the addition of DICOM Query/Retrieve.

StealthStation® S7® System Specifications

- Surgeon Monitor cart footprint (at widest point): 23.75" x 24.86" (60.3 cm x 63.1 cm)
- Staff Cart footprint (at widest point): 23.01" x 23.84" (58.4 cm x 60.5 cm)
- Surgeon Monitor viewable size: 24" widescreen (61 cm)
- Surgeon Monitor resolution: 1920 x 1200, 60 Hz
- Surgeon Monitor option: Equipped with or without touchscreen interface
- Dedicated staff monitor and workspace
- Staff Monitor viewable size: 15.4" widescreen (39.1 cm)
- Staff Monitor resolution: 1440 x 900, 60 Hz
- Maximum camera height: 94.5" (240 cm)
- Maximum camera reach: 49.2" (125 cm)
- AxiEM™ compatible – industry's only electromagnetic option for neuronavigation
- O-arm® System, Polestar® System, and intra-operative imaging system compatible
- Microscope interface integration with most manufacturers
- DICOM Query/Retrieve
- USB Stick or CD/DVD transfer of exams to the StealthStation® S7® System
- High Fidelity MP3 Audio Speaker System (MP3 player not included)

IMPROVED WORKFLOW

Operating Rooms are dynamic work environments where everyone plays a key role in providing care for the patient. That is why every aspect of the StealthStation® S7® was thoughtfully considered. You deserve a friendly experience throughout set up, navigation, shut down, and transport.

- **Intelligent, 2-Cart Design.** The Surgeon monitor cart undocks with one hand and has an agile spring arm for optimal positioning of the wide screen monitor around the surgical team.
- **Compact Footprint.** OR floor space is critical real estate. The StealthStation® S7® AxiEM™ system's sleek, upright form compresses footprint without compromising functionality.
- **MP3 Music Ready.** We've equipped the StealthStation® S7® System with a high fidelity sound system to make your OR environment just a little more enjoyable.
- **Attention to Details.** Meticulous consideration of every detail helps make the StealthStation S7® System an intuitive and advanced navigation system.
 - The StealthStation® S7® System works in the docked or undocked positions.
 - The AxiEM™ tracking system can be used from its discreet location on a StealthStation® S7® System or easily lifted and attached to the OR bedrail.
 - A longer camera handle and laser button make it is easier to aim the navigation camera from any position.
 - The StealthStation® S7® system includes a roomy storage drawer, equipment bin, and cable wraps to keep things organized without taking up more space.
 - The latest in uninterrupted power systems allows you to unplug and move the StealthStation® S7® System without losing data.

www.medtronicneurosurgery.com

Medtronic, Inc.

Surgical Technologies, Neurosurgery
826 Coal Creek Circle, Louisville, CO 80027 USA
Toll Free: 877-242-9504 | Telephone: 720-890-3200 | Fax: 720-890-3500



Rx only. Refer to product instruction manual/package insert for instructions, warning, precautions and contraindications. For further information, please contact Medtronic ST Neurosurgery at (877) 242-9504, and/or consult Medtronic's website at www.medtronicneurosurgery.com.

©2013 Medtronic, Inc. All Rights reserved. Printed in the USA. StealthStation®, S7®, O-arm® and PoleStar® are registered trademarks of Medtronic, Inc. AxiEM™ is a trademark of Medtronic, Inc. UC201402388 EN 11.13